


# ezject

electronic injector

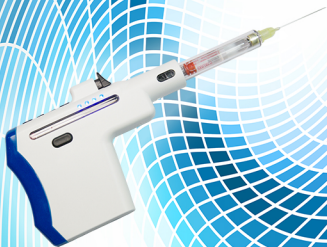

## **Painless Injection**

Fine speed control of injection minimize pain

## **Periodontal Injection**

Designed to inject the needle easily

## **Infiltration Anesthesia**

Local anesthesia is possible

Wellness of Your Life!

**DIGIMED**

#### 4 step Control of Injection Speed

As Dental syringe, ezject is designed to minimize patient's pain by keeping stable syringe speed.

#### Ampule Indicator

Injection amount is easily identified.

#### 3 step Pressure Sensor

Additional speed control with a pressure sensor



#### Injection Speed Indicating LED

LED indicates the injection speed by 4 steps.

#### Battery Recharging Dock

Using rechargeable battery, ezject is easy to move and bring anywhere.

### Specification

Model Name	ezject
Injection Control Method	Micro-processor and Software control
Maximum Injection Distance	50 mm
Injection Speed Step	4 step
Injection Speed Control	6 step control by magnetic force
Function	Aspiration features Injection amount check with analog indicator Battery level check Auto power-off LED injection speed indicator Auto location after full injection
Total Weight	130 g (ezject)
Dimension	21(w) x 123(d) x 88(h)
Power Requirement	Battery: 7.4 V / Charger: AC 100 ~ 240 V

Distributed by

#### DIGIMED CO., LTD

309~311-ho  
145, Gasan digital 1-ro, Geumcheon-gu  
Seoul 08506 Rep. of Korea  
TEL 82-2-2624-1551  
FAX 82-2-2624-1553  
digimed@digimed.co.kr  
www.digimed.co.kr